

# Femur Nailing System

For Secure and Stable fixation.

**Miraclus**<sup>®</sup>  
Implanting LIFE



# Femur Nailing System

For Secure and Stable fixation.

## Femur Nail Indication:


- Indicated for Femoral shaft fracture.

## Features:

**A. EntryPoint:** Tip of the greater trochanter.


## B. Proximal Locking:

- I. **Static Locking** - Static locking is done with 5 mm locking bolt at 120° antegrade position.

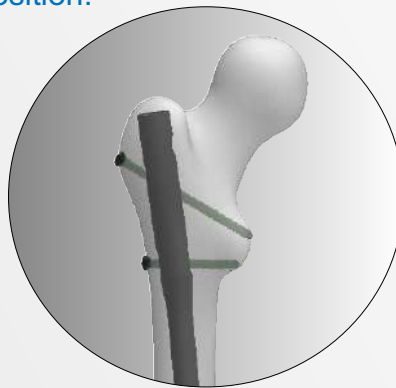


- II. **Dynamic Locking** -

- a. The dynamic locking is possible at proximal position of the dynamic capsule(hole) slot.

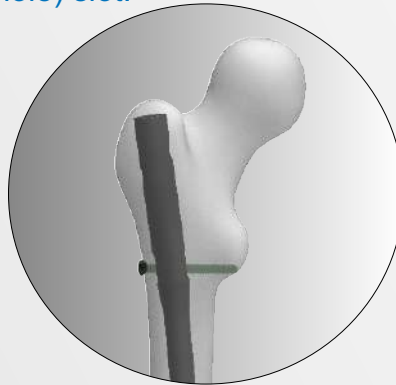

- b. The static locking is possible at distal position of the standard slot.

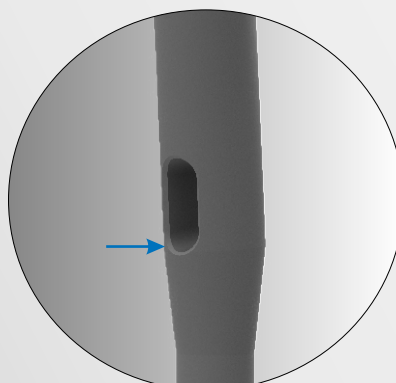


# Femur Nailing System


For Secure and Stable fixation.

## B. Distal Locking:

I. **Static Locking** - static locking is done with 5 mm locking bolt in hole A & B.

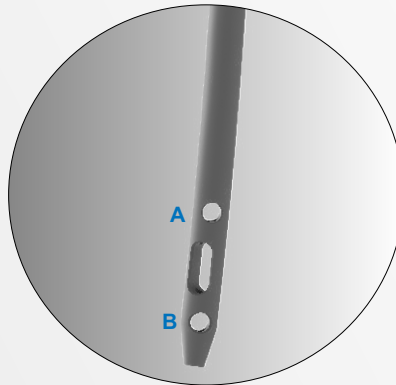

II. **Dynamic Locking** - dynamic locking is done with 5 mm at distal position of the dynamic capsule hole.

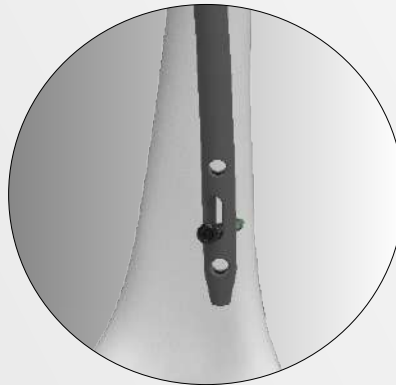

## Available Femur Nail Sizes :

Dia (mm)	Length (mm)				
	340	360	380	400	420
9 (L & R)	340	360	380	400	420
10 (L & R)	340	360	380	400	420
11 (L & R)	340	360	380	400	420